


## Zone 70 - Contaminant Lab Allergens and Microbiology

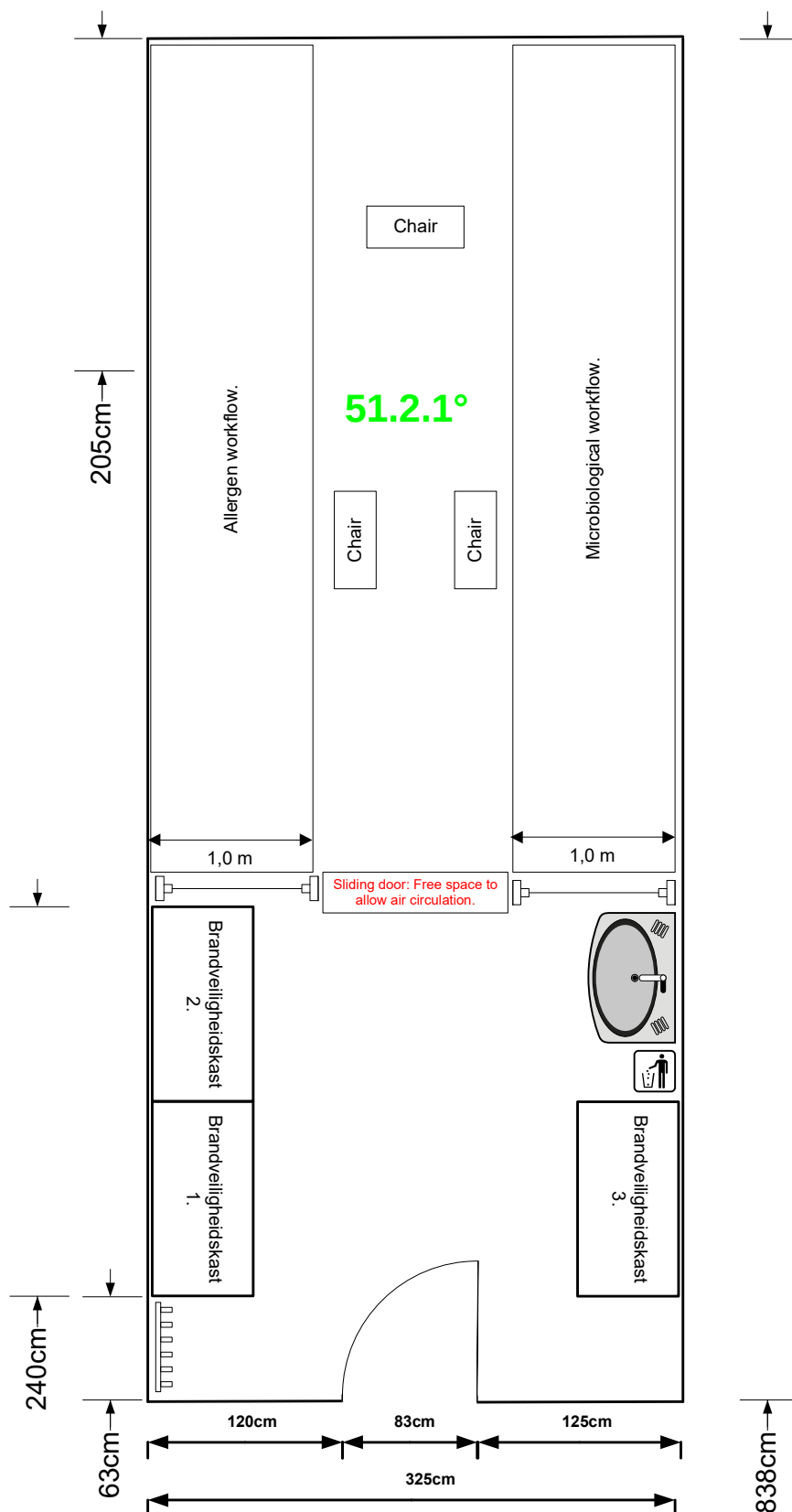

### Exact dimensions:

Width area: 3m25.  
Depth area: 8m38.  
Door: 83 cm.  
Door left: 119 cm.  
Door right: 123 cm.  
Fire safety cabinet: 62.5 cm \* 119.5 cm.

Available space: 5m35 \* 3m25 (17.4 m2).

### Required utilities and modifications:

**Sink:** Cold and hot water beside table to avoid moisture contamination.

**Laminar Flow:** Air introduced within area must pass laminar flow.

Required maintenance laminar flow to be requested.

**Storage cabinets:** Left and right of laminar flow.

**Tables:** Equipped drawers and cabinets.  
**Refrigerator.**

**Power connections:** See below.

**Entrance door:** To be modified with window or replaced by transparent door.

Separation entrance - analytical part recommended.

**UV-light:** Activation when outside door locked.

**Focus point:** Waste removal. Area must be locked when not in operation.

### Key equipment Allergen workflow requiring power: 5 PC - 2 LC.

Reader.  
Elisa equipment.

### Key equipment Microbiology workflow requiring power: 15 PC - 2 LC.

Laminar flow.  
Stomacher ("Paddle" blender).  
Incubators (20°C / 37°C / 41°C).  
Autoclave.  
Dry bath (2).  
MD system.  
PC / Screen.  
UV lamps: 2.

Lan connections (LC): 4.  
Power connections (PC): 20.

## Zone 70 - Contaminant Lab Allergens and Microbiology - Details

